



Developed by The Specialist in Filtration

The Ultimate YAMASHIN FILTER MASK

フィルタ専門メーカーが開発した


究極の ヤマシン・ フィルタマスク



Bring a medical grade mask to your home! "Protect your loved ones"

Combine nano-level filters featuring three-dimensional-multi-layer structure, which are much finer than the micro levels, to protect you and your people from invisible threats.

- World-class Filtration Technology
- Innovative 3D Nano Technology Filter
- 100% Made in Japan
- Minimal deterioration of filtration efficiency over 24 hours

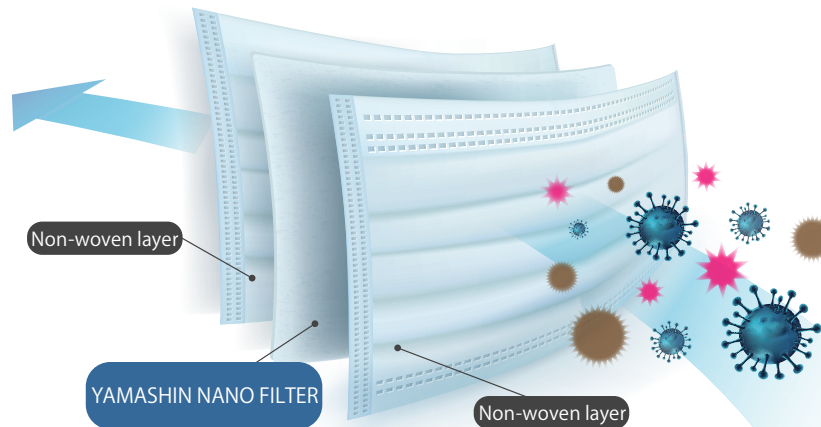


Unlike disposable masks that trap particles by electrostatic charge, YAMASHIN FILTER MASK catches particles with the fineness of the filter itself, which is 3D structured by our patented special nano fiber technology. Thus, it does not deteriorate as quickly over 24 hours.

*YAMASHIN FILTER MASK has filtration efficiency equivalent to N95, a standard by U.S. National Institute for Occupational Safety and Health, on the ground that the filter media are shared by another mask model we have for medical use, which is in process for DS2 certification by Japanese Government, said to be equivalent to N95.

Block airborne droplet by the nano level multi-layer


Revolutionary filtration technology independent of electrostatic charge
super fine nano fibers

YAMASHIN FILTER MASK
Fiber diameter: 200nm(0.2μm)

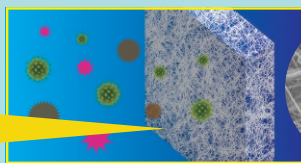

Electronic microscope X1000

Generic mask filter
Fiber diameter Avg. 3μm

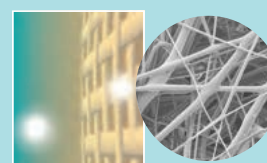


Electronic microscope X1000


Generic masks depend on electrostatic charge for filtration efficiency. Electrostatic charge deteriorates when exposed to air or humidity

99% of 0.1μm Fine Particles Captured by the YAMASHIN NANO FILTER

※Results of PFE tests conducted by Kaken Test Center (Japanese third party test organization)

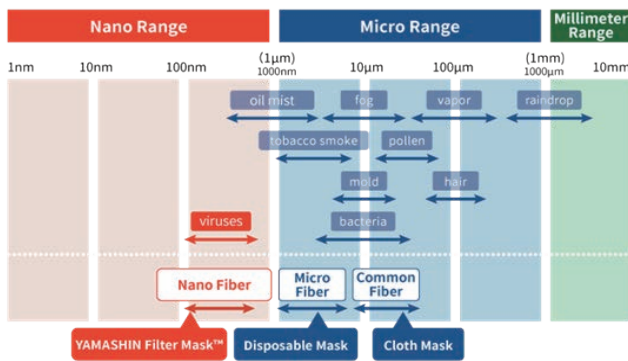

High-performance filter made of nano fibers keeps the initial efficiency even after 5-time hand washes

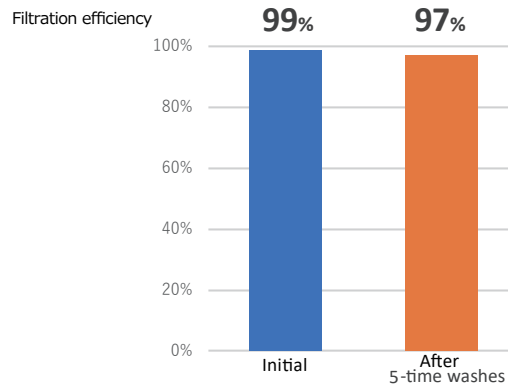

※Avg. of PFE test results conducted by Kaken Test Center
※Washed by hand

